

Wash-Check cleaning efficacy system for Ultrasonics

Valisafe
clearly a better indicator

Wash-Checks U Ultrasonic Monitor

LOT 090000

0000-00



Easy to use



- For monitoring of the cleaning efficacy of Ultrasonic cleaning systems
- Test soil is designed to mimic the removal of blood and tissue during the Ultrasonic process.
- Low cost to ensure routine use
- Easy to read
- Stainless steel holder (for use with both the Valisafe standard wash checks and the Valisafe Wash check U)
- Indicators are non toxic and have a 2 year shelf life
- Durable stainless steel holder

Valisafe Wash-Check cleaning efficacy system for Ultrasonics


The Valisafe Wash-check indicators for Ultrasonics when used in conjunction with the wash check holder is designed to establish an effective routine monitoring solution to monitor the ultrasonic soil removal efficacy. The Valisafe Wash-check U is recommended at the beginning of every day in an empty load to provide machine release and then periodically throughout the day with instruments for load release.

Order Details


V3503950
V3503930

Wash-Checks U
Wash-Check holder


50 indicator / bag – 5 bags / case
6 holders per case

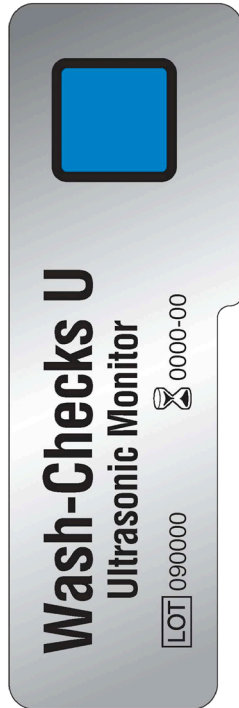


Name:	Wash-Checks ultrasonic cleaning efficacy Wash-Checks holder
Classification:	AAMI ISO 15883
Storage conditions:	room temperature in sealed foil bag
Lot no:	yes
Prod date:	yes
Exp date:	yes
Correct symbols:	In accordance with AAMI ISO 15883
Indicator:	Cleaning efficacy indicator
Non toxic:	yes
Disp/single use::	indicator single use - holder multiple uses

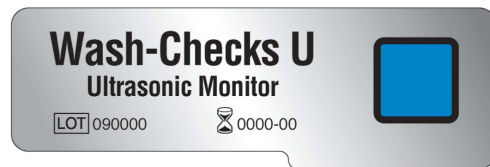

Exposed

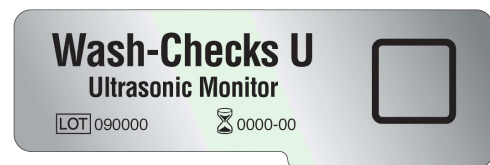

Correct colour change

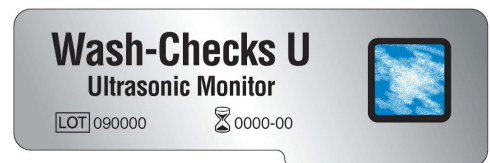

Insufficient